

# CROOKES TUBE

N99-B10-7260

\$695

Used for experiments and observations of the deflection of cathode rays by an electric or magnetic field. Can also be used to explain the principle of a cathode ray tube (CRT). The incorporation of a hot-cathode, gas filled electron tube provides a clear view of the electron beam.




## SPECIFICATIONS

Tube: 75mm, equipped with terminal block for connection to power supply and support stand.

Power supply: Heater: 6.3V/2A Cathode Plate: 0-400V DC

Dimensions: 305x95x110mm

Weight: 350g

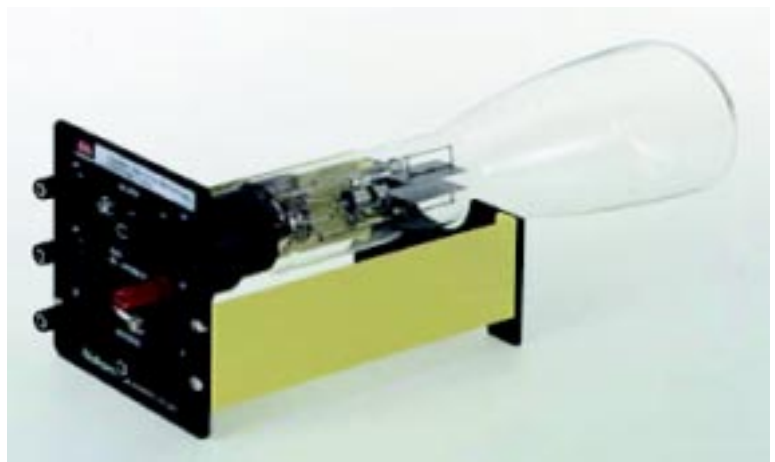

---

# POWER SUPPLY FOR CROOKES TUBE

N50-1000

\$389



## SPECIFICATIONS

- Tube: 75mm, equipped with terminal block for connection to power supply and support stand
- Power Supply: Heater: 6.3V/2A, Cathode plate: 0-400V DC, Dimensions: 305x95x110mm
- Weight: 350g

This power supply unit is specially designed for use with the Crookes Tube.

Call 1-800-799-6232, shop online at [nadascientific.com](http://nadascientific.com) or contact your favorite science dealer today.