
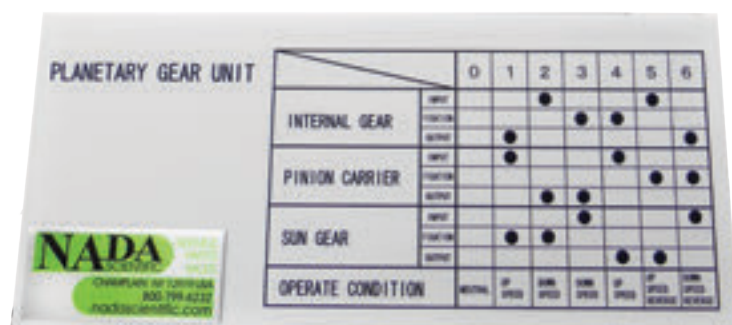


## PLANETARY GEAR MODEL

N98-V061-SPG

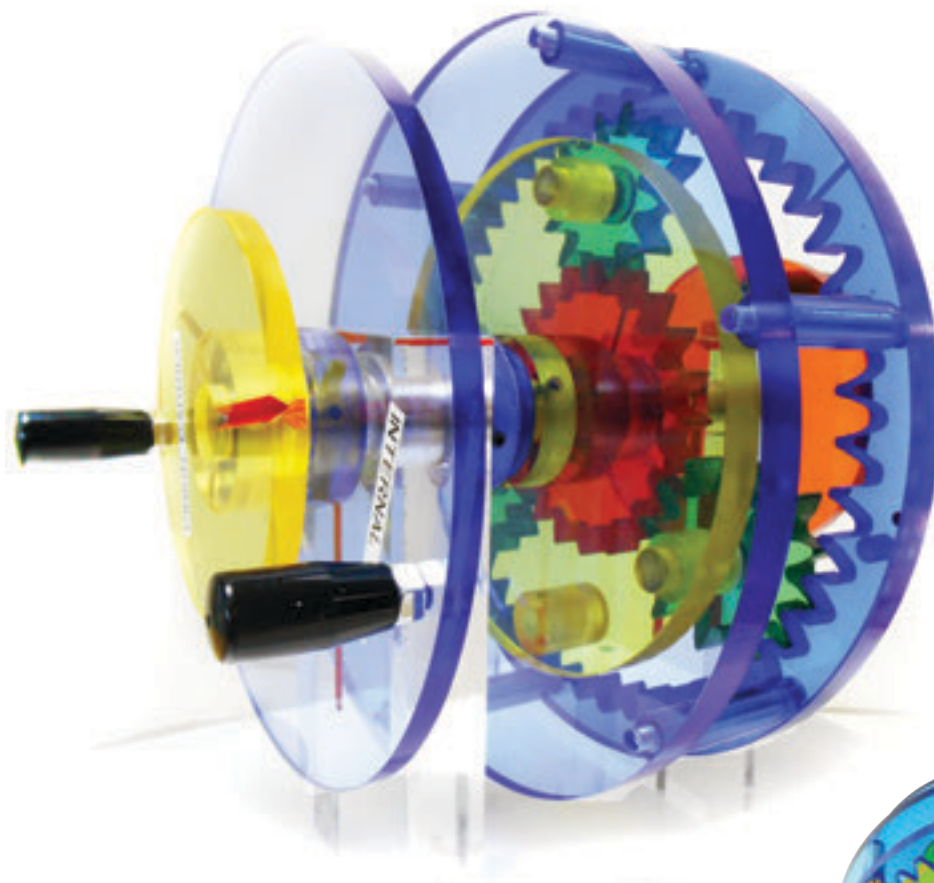

This technical model acts similarly to the Ravigneaux planetary gearset of a modern vehicle. The complex principles of a planetary gear are simplified by the use of color coded discs. Blue is used for the internal gear. Green is used for the pinion gear. Yellow is used for the sun gear. This model displays 3 forward speeds and 1 reverse speed. You may operate this unit by locking down a specific part of the gear using a fixed pin. Changes in rotation are visible when adjusting either the sun, pinion, or internal gear. A nomenclature panel is included on the front of the model for reference.

### Weight & Dimensions

27cm x 30cm x 24cm (L x W x H)

Net Weight: 10 kg

MOUNTED ON  
BASE

Version2.2SB091119

